


Qymia Workflow Example

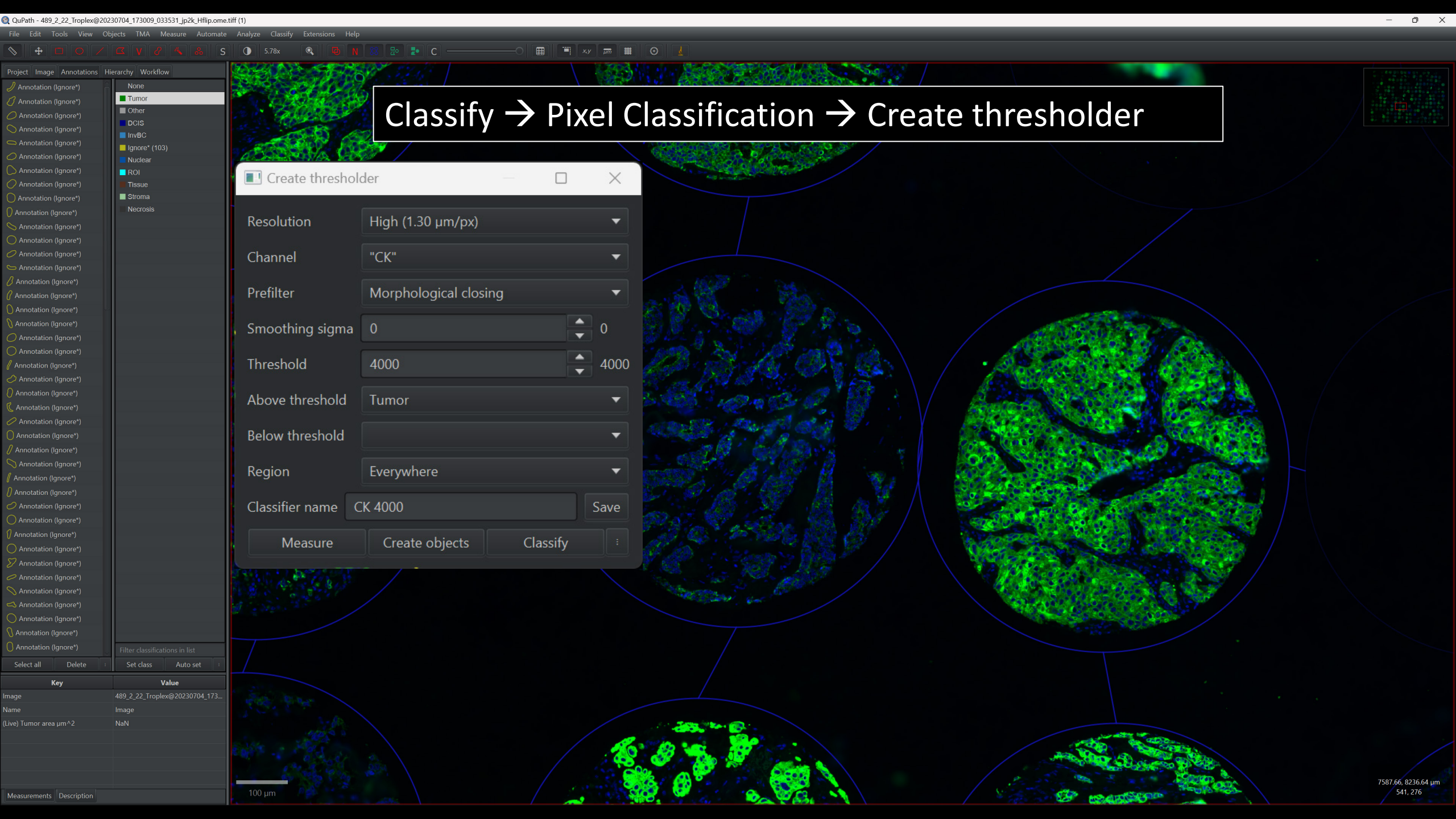

Classify → Pixel Classification → Create threshold

Create threshold

Resolution

High (1.30 μm/px)

Channel

"CK"

Prefilter

Morphological closing

Smoothing sigma

0

0

Threshold

4000

4000

Above threshold

Tumor

Below threshold

Region

Everywhere

Classifier name

CK 4000


Save

Measure

Create objects

Classify

Key	Value
Image	489_2_22_Troplex@20230704_173...
Name	Image
(Live) Tumor area μm^2	NaN
Measurements	Description

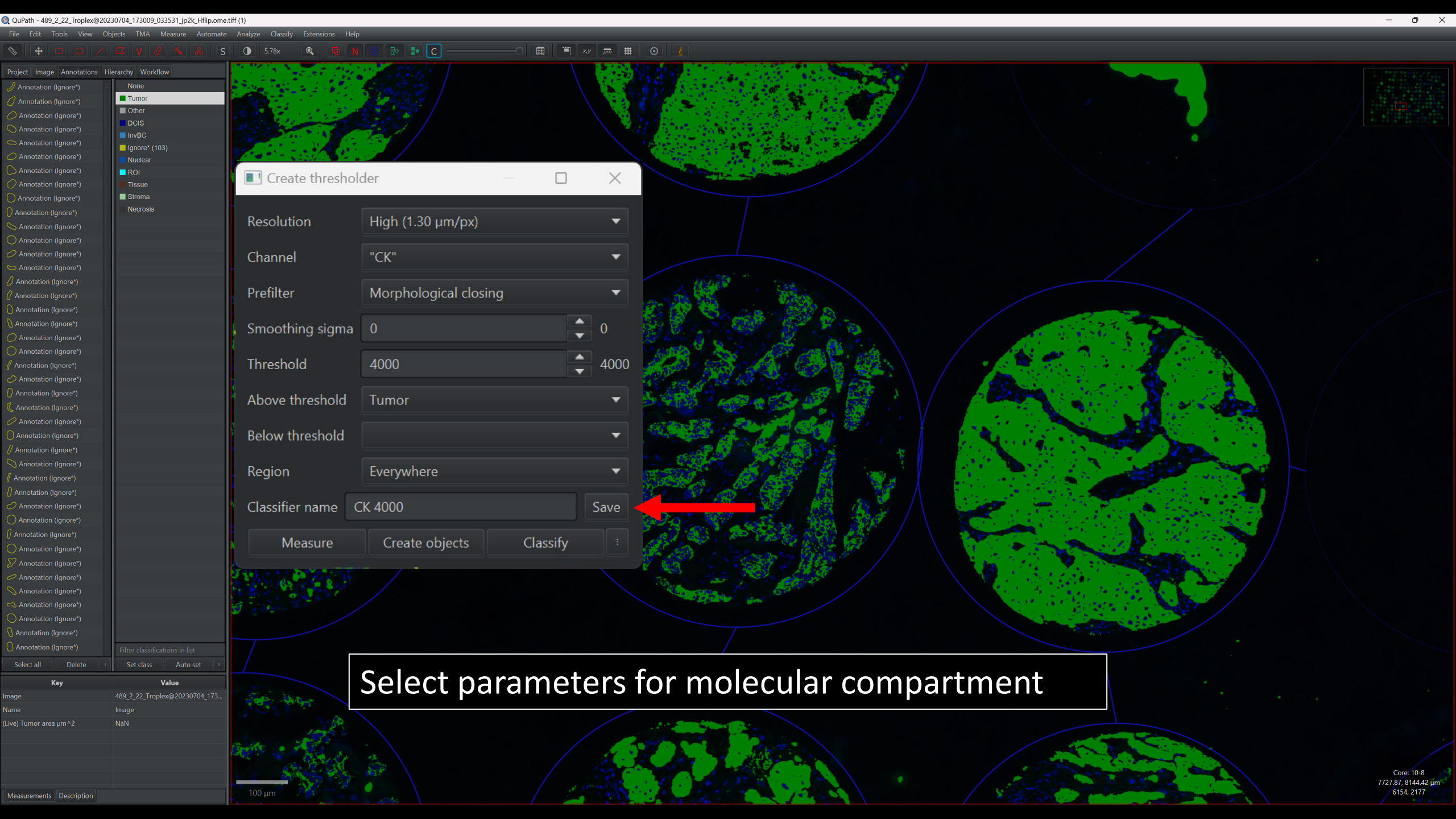

Create threshold

Resolution

High (1.30 µm/px)

Channel

"CK"

Prefilter

Morphological closing

Smoothing sigma

0

0

Threshold

4000

4000

Above threshold

Tumor

Below threshold

Region

Everywhere

Classifier name

CK 4000

Save

Measure

Create objects

Classify

Select parameters for molecular compartment

Key	Value
Image	489_2_22_Troplex@20230704_173...
Name	Image
(Live) Tumor area µm^2	NaN
Measurements	Description

Project Image Annotations Hierarchy Workflow

Create objects

?

New object type: Annotation

Minimum object size: 5 μm^2

Minimum hole size: 0 μm^2

☐ Split objects

☒ Delete existing objects

☐ Create objects for ignored classes

☐ Set new objects to selected

OK Cancel

Core: 10-8
7727.87, 8144.42 μm
6154, 2177


The screenshot displays the QuPath software interface for digital pathology analysis. The main workspace shows a series of histological images with segmentation masks applied. A legend on the left side lists the following classification classes:

- Tumor (186)
- Other
- DCIS
- InvBC
- Ignore* (103)
- Nuclear
- ROI
- Tissue
- Stroma
- Necrosis

A "Check objects" dialog box is currently open, allowing users to verify or modify the segmentation results. The bottom left corner features a table with metadata for the loaded image:

Key	Value
Image	489_2_22_Tropex@20230704_173...
Name	Image


The bottom right corner provides additional information, including the core ID (Core: 10-8), coordinates (776855, 8261.64 μm), and pixel dimensions (6996, 1633). A scale bar indicating 100 μm is visible in the bottom left corner of the image area.

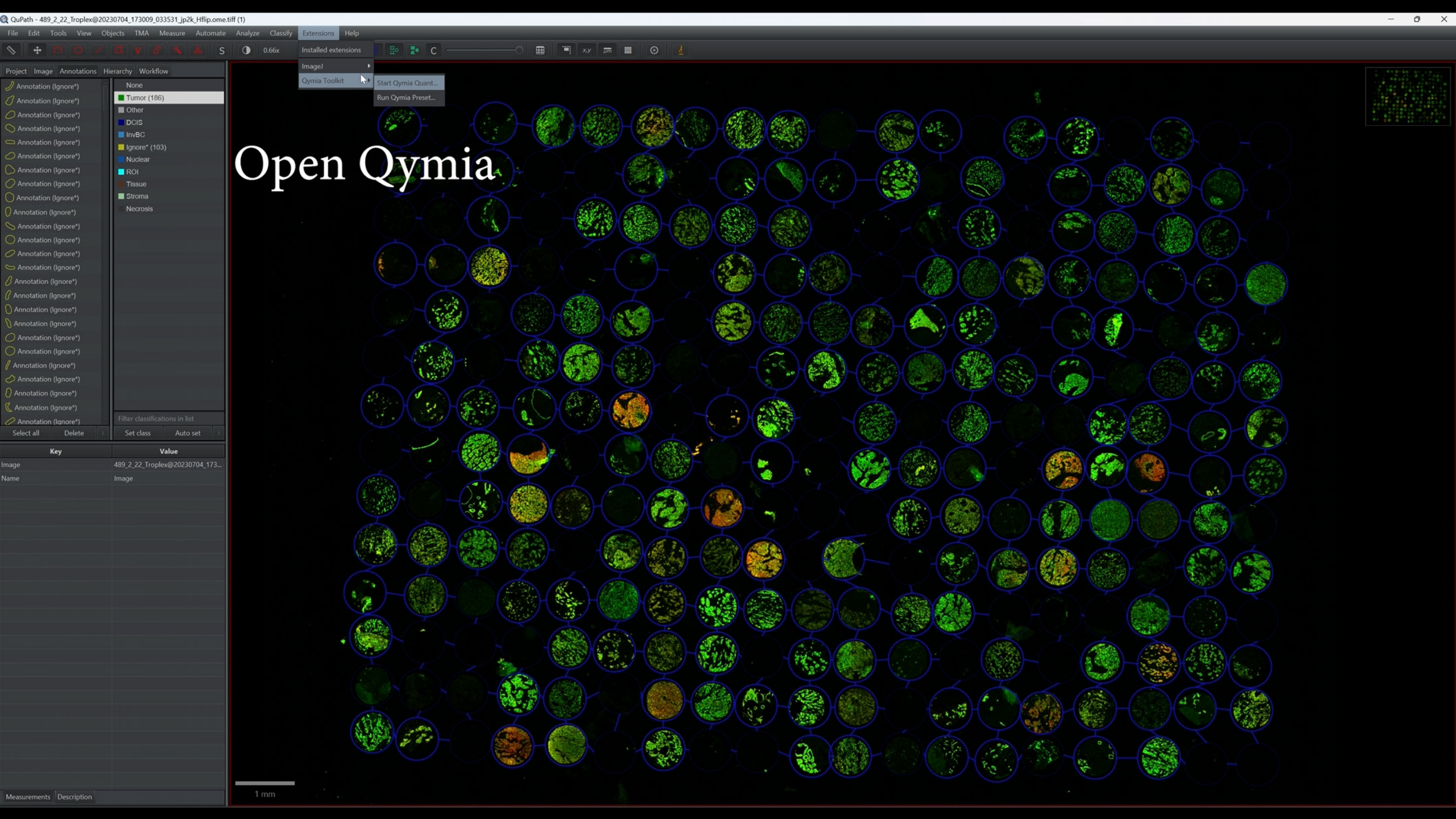


Open Qymia

Key	Value
Image	489_2_22_Troplex@20230704_173...
Name	Image

Set imaging parameters

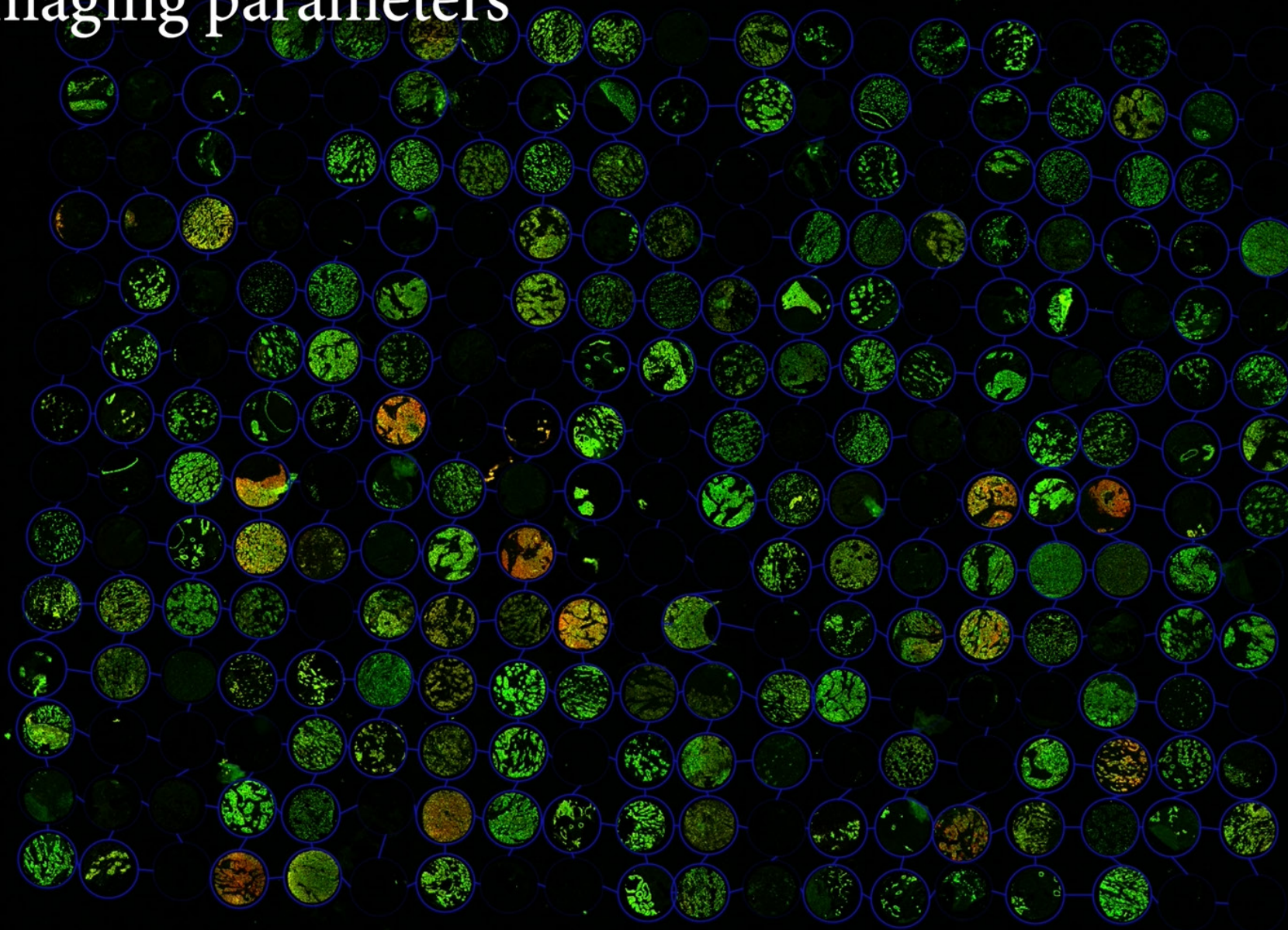

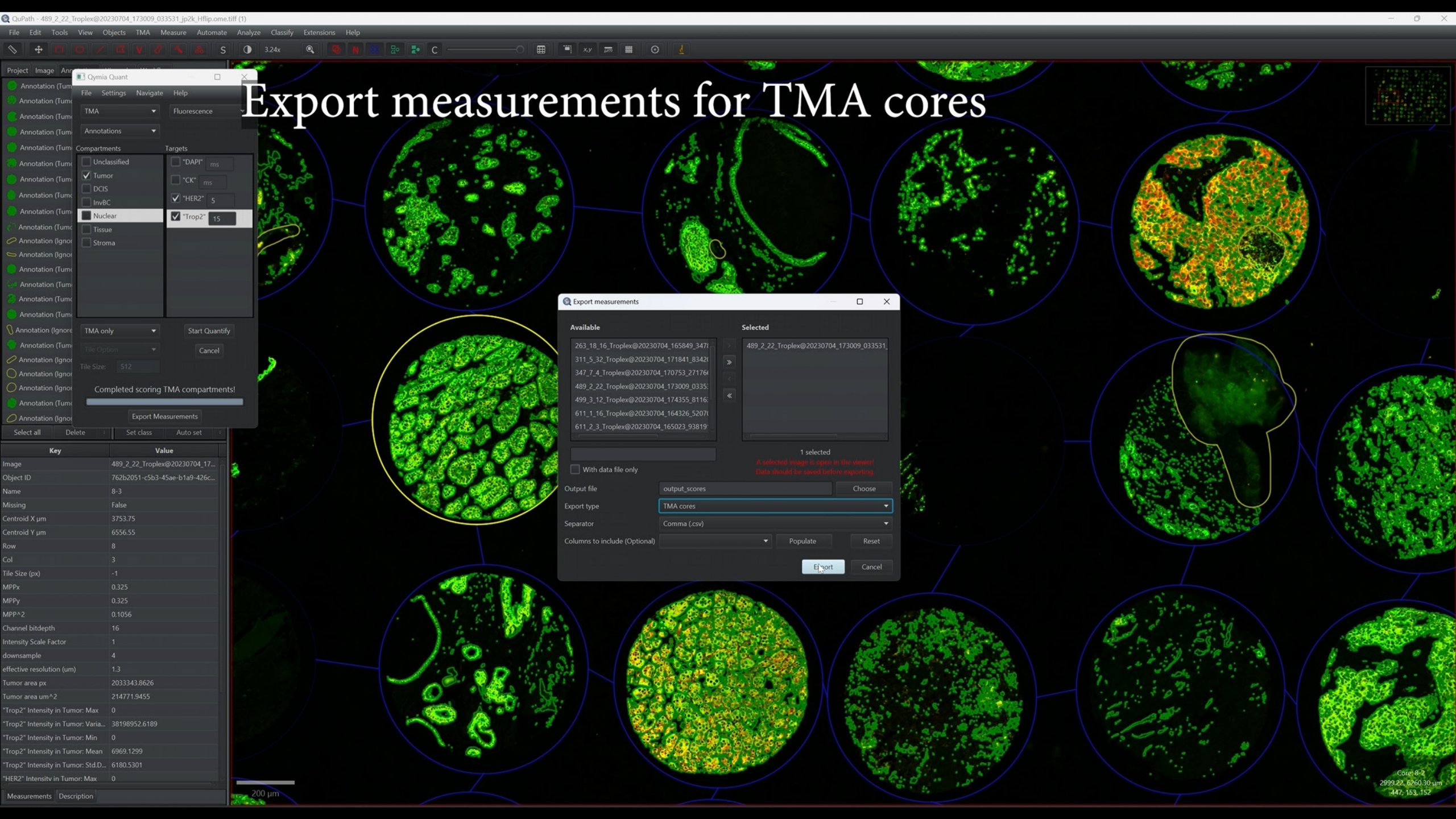

Export measurements for TMA cores

Key	Value
Image	489_2_22_Troplex@20230704_17...
Object ID	762b2051-c5b3-45ae-b1a9-426c...
Name	8-3
Missing	False
Centroid X μm	3753.75
Centroid Y μm	6556.55
Row	8
Col	3
Tile Size (px)	-1
MPPx	0.325
MPPy	0.325
MPP $\wedge$ 2	0.1056
Channel bitdepth	16
Intensity Scale Factor	1
downsample	4
effective resolution (μm)	1.3
Tumor area px	2033343.8626
Tumor area μm^2	214771.9455
"Trop2" Intensity in Tumor: Max	0
"Trop2" Intensity in Tumor: Varia...	38198952.6189
"Trop2" Intensity in Tumor: Min	0
"Trop2" Intensity in Tumor: Mean	6969.1299
"Trop2" Intensity in Tumor: Std.D...	6180.5301
"HER2" Intensity in Tumor: Max	0
Measurements	Description

Core: 8-3
2999.21, 6260.30 μm
447, 153, 152